



(12) **United States Patent**  
**Behzadi**

(10) **Patent No.:** **US 10,463,505 B2**  
(45) **Date of Patent:** **Nov. 5, 2019**

(54) **BONE PREPARATION APPARATUS AND METHOD**

(56) **References Cited**

U.S. PATENT DOCUMENTS

(71) Applicant: **Kambiz Behzadi**, Pleasanton, CA (US)  
(72) Inventor: **Kambiz Behzadi**, Pleasanton, CA (US)  
(\*) Notice: Subject to any disclaimer, the term of this patent is extended or adjusted under 35 U.S.C. 154(b) by 255 days.

3,818,514 A	6/1974	Clark
4,530,114 A	7/1985	Tepic
4,712,951 A	12/1987	Brown
5,133,765 A	7/1992	Cuilleron
5,318,570 A	6/1994	Hood et al.
5,534,006 A	7/1996	Szabo et al.
5,769,092 A	6/1998	Williamson, Jr.
5,849,015 A	12/1998	Haywood et al.
6,110,179 A	8/2000	Flivik et al.
6,161,545 A *	12/2000	Chow ..... A61B 17/320068 128/898
6,585,771 B1	7/2003	Buttermilch et al.
7,036,211 B1	5/2006	Panks

(Continued)

(21) Appl. No.: **15/235,053**

(22) Filed: **Aug. 11, 2016**

(65) **Prior Publication Data**

US 2017/0196705 A1 Jul. 13, 2017

**Related U.S. Application Data**

(63) Continuation-in-part of application No. 15/202,434, filed on Jul. 5, 2016.  
(60) Provisional application No. 62/277,294, filed on Jan. 11, 2016.

(51) **Int. Cl.**

**A61F 2/46** (2006.01)  
**A61B 17/16** (2006.01)  
**A61B 17/32** (2006.01)  
**A61B 17/14** (2006.01)

(52) **U.S. Cl.**

CPC ..... **A61F 2/4603** (2013.01); **A61B 17/142** (2016.11); **A61B 17/1666** (2013.01); **A61B 17/320068** (2013.01); **A61F 2/4607** (2013.01); **A61F 2/4609** (2013.01); **A61F 2/4612** (2013.01); **A61F 2002/4627** (2013.01); **A61F 2002/4681** (2013.01); **A61F 2002/4683** (2013.01)

(58) **Field of Classification Search**

CPC ..... A61F 2/46; A61F 2/4603; A61F 2/4609; A61F 2/4607; A61F 2/4612; A61B 17/16; A61B 17/1666

See application file for complete search history.

FOREIGN PATENT DOCUMENTS

EP	1433445 A1	6/2004
WO	2018031752 A1	2/2018

OTHER PUBLICATIONS

International Search Report for International application No. PCT/US2017/046261, dated Oct. 18, 2017.

(Continued)


*Primary Examiner* — Christopher J Beccia

(74) *Attorney, Agent, or Firm* — Patent Law Offices of Michael E. Woods; Michael E. Woods

(57) **ABSTRACT**

A system and method for improving installation of a prosthesis. Devices include prosthesis installation tools, prosthesis assembly tools, site preparation systems, and improved power tools used in implant site preparation.

**15 Claims, 15 Drawing Sheets**

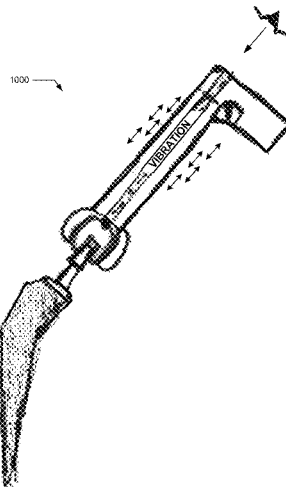

(56)

## References Cited

## U.S. PATENT DOCUMENTS

8,167,823	B2	5/2012	Nycz et al.	
8,603,100	B2	12/2013	Muller	
2002/0183851	A1	12/2002	Spiegelberg et al.	
2003/0065398	A1	4/2003	Cueille et al.	
2003/0229357	A1	12/2003	Dye	
2005/0101962	A1	5/2005	Schwenke et al.	
2006/0015110	A1*	1/2006	Pepper	A61B 17/1615 606/80
2006/0142754	A1	6/2006	Irion et al.	
2007/0162038	A1	7/2007	Tuke	
2009/0192626	A1	7/2009	Keefter et al.	
2010/0249796	A1	9/2010	Nycz	
2011/0004318	A1	1/2011	Tulkis et al.	
2011/0178521	A1	7/2011	Siravo et al.	
2011/0264009	A1	10/2011	Walter et al.	
2012/0172939	A1*	7/2012	Pedicini	A61B 17/1604 606/86 R
2012/0209277	A1	8/2012	Leparmentier et al.	
2013/0211535	A1	8/2013	Cueille	
2013/0226189	A1	8/2013	Young	
2014/0135773	A1	5/2014	Stein et al.	
2014/0135791	A1	5/2014	Nikou et al.	
2015/0182350	A1	7/2015	Behzadi	
2015/0182351	A1	7/2015	Behzadi	
2015/0196343	A1	7/2015	Donald et al.	
2015/0201918	A1	7/2015	Kumar et al.	
2015/0282856	A1	10/2015	Haiat et al.	
2016/0166390	A1	6/2016	Dye et al.	
2016/0206430	A1	7/2016	Grostefon et al.	
2016/0206433	A1	7/2016	Grostefon et al.	
2016/0220315	A1	8/2016	Falardeau et al.	
2017/0056205	A1	3/2017	Biegun et al.	
2017/0196506	A1	7/2017	Behzadi	
2017/0196701	A1	7/2017	Behzadi et al.	
2017/0196704	A1	7/2017	Behzadi et al.	
2017/0196705	A1	7/2017	Behzadi	
2017/0196706	A1	7/2017	Behzadi	
2017/0196707	A1	7/2017	Behzadi	
2017/0196708	A1	7/2017	Behzadi et al.	
2017/0196710	A1	7/2017	Behzadi	
2017/0196711	A1	7/2017	Behzadi	
2017/0325972	A1	11/2017	Steif	
2017/0340456	A1	11/2017	Behzadi	
2018/0235764	A1	8/2018	Moore et al.	

## OTHER PUBLICATIONS

Written Opinion of the International Searching Authority for International application No. PCT/US2017/046261, dated Oct. 18, 2017.

International Search Report for International application No. PCT/US2017/012753, dated May 5, 2017.

Written Opinion of the International Searching Authority for International application No. PCT/US2017/012753 dated May 5, 2017.

U.S. Appl. No. 16/276,639, filed Feb. 15, 2019, Kambiz Behzadi.

U.S. Appl. No. 16/278,085, filed Feb. 16, 2019, Kambiz Behzadi.

U.S. Appl. No. 16/278,668, filed Feb. 18, 2019, Kambiz Behzadi.

U.S. Appl. No. 16/374,750, filed Apr. 4, 2019, Kambiz Behzadi et al.

U.S. Appl. No. 16/375,736, filed Apr. 4, 2019, Kambiz Behzadi et al.

U.S. Appl. No. 62/277,294, filed Jan. 11, 2016, Kambiz Behzadi.

U.S. Appl. No. 62/353,024, filed Jun. 21, 2016, Kambiz Behzadi.

U.S. Appl. No. 62/355,657, filed Jun. 28, 2016, Kambiz Behzadi.

U.S. Appl. No. 62/373,515, filed Aug. 11, 2016, Kambiz Behzadi.

U.S. Appl. No. 62/651,077, filed Mar. 31, 2018, Kambiz Behzadi.

U.S. Appl. No. 62/742,851, filed Oct. 8, 2018, Kambiz Behzadi.

U.S. Appl. No. 62/743,042, filed Oct. 9, 2018 Kambiz Behzadi et al.

U.S. Appl. No. 15/202,434, filed Jul. 5, 2016, Kambiz Behzadi.

U.S. Appl. No. 15/234,782, filed Aug. 11, 2016, Kambiz Behzadi et al.

U.S. Appl. No. 15/234,880, filed Aug. 11, 2016, Kambiz Behzadi et al.

U.S. Appl. No. 15/235,032, filed Aug. 11, 2016, Kambiz Behzadi et al.

U.S. Appl. No. 15/284,091, filed Oct. 3, 2016, Kambiz Behzadi.

U.S. Appl. No. 15/362,675, filed Nov. 28, 2016, Kambiz Behzadi.

U.S. Appl. No. 15/396,785, filed Jan. 2, 2017, Kambiz Behzadi et al.

U.S. Appl. No. 15/398,996, filed Jan. 5, 2017, Kambiz Behzadi.

U.S. Appl. No. 15/453,219, filed Mar. 8, 2017, Kambiz Behzadi.

U.S. Appl. No. 15/592,229, filed May 11, 2017, Kambiz Behzadi.

U.S. Appl. No. 15/687,324, filed Aug. 25, 2017, Kambiz Behzadi.

U.S. Appl. No. 15/716,529, filed Sep. 27, 2017, Kambiz Behzadi et al.

U.S. Appl. No. 15/716,533, filed Sep. 27, 2017, Kambiz Behzadi.

U.S. Appl. No. 16/030,603, filed Jul. 9, 2018, Kambiz Behzadi.

U.S. Appl. No. 16/030,824, filed Jul. 9, 2018, Kambiz Behzadi.

U.S. Appl. No. 16/154,033, filed Oct. 8, 2018, Kambiz Behzadi et al.

\* cited by examiner


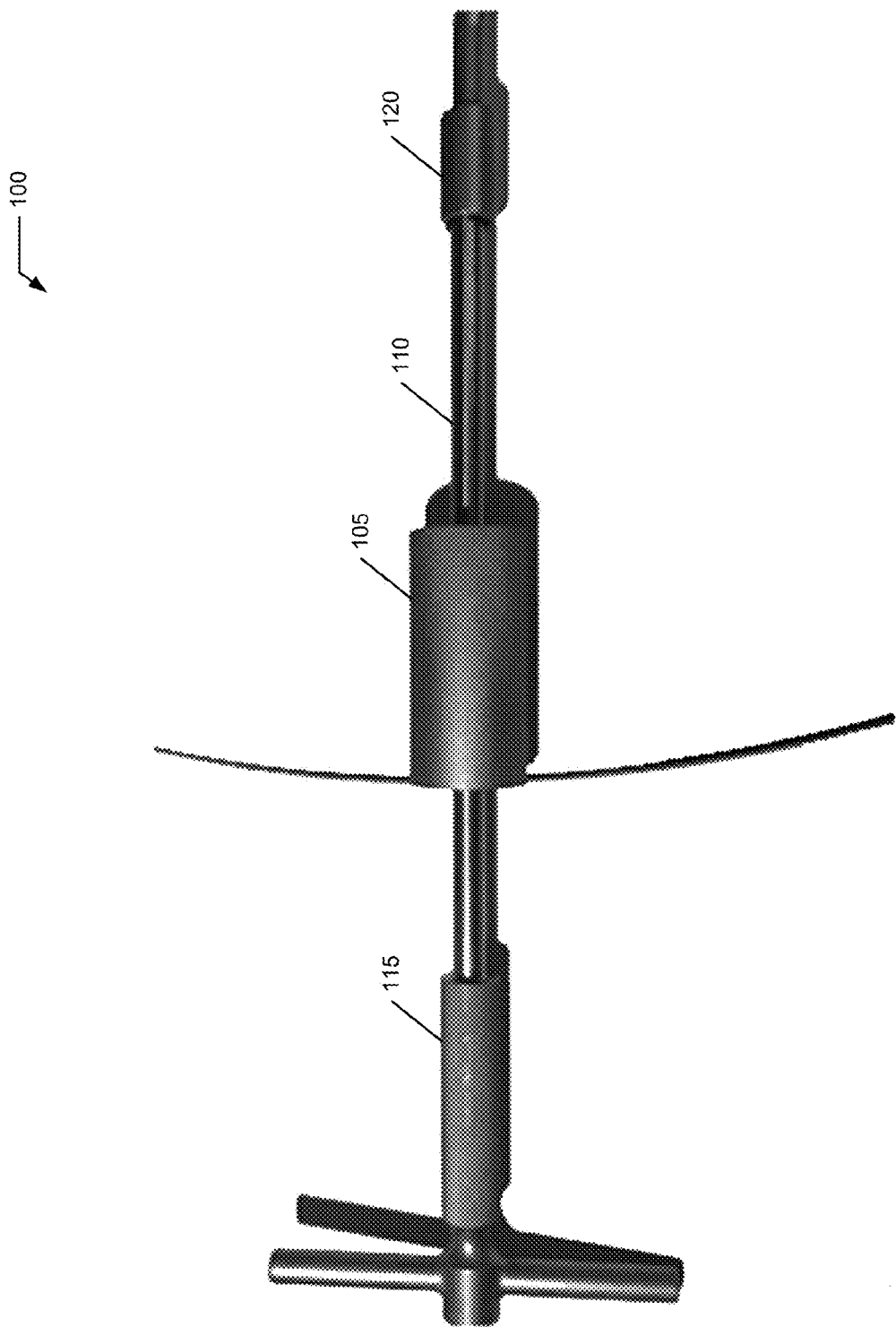


FIG. 1

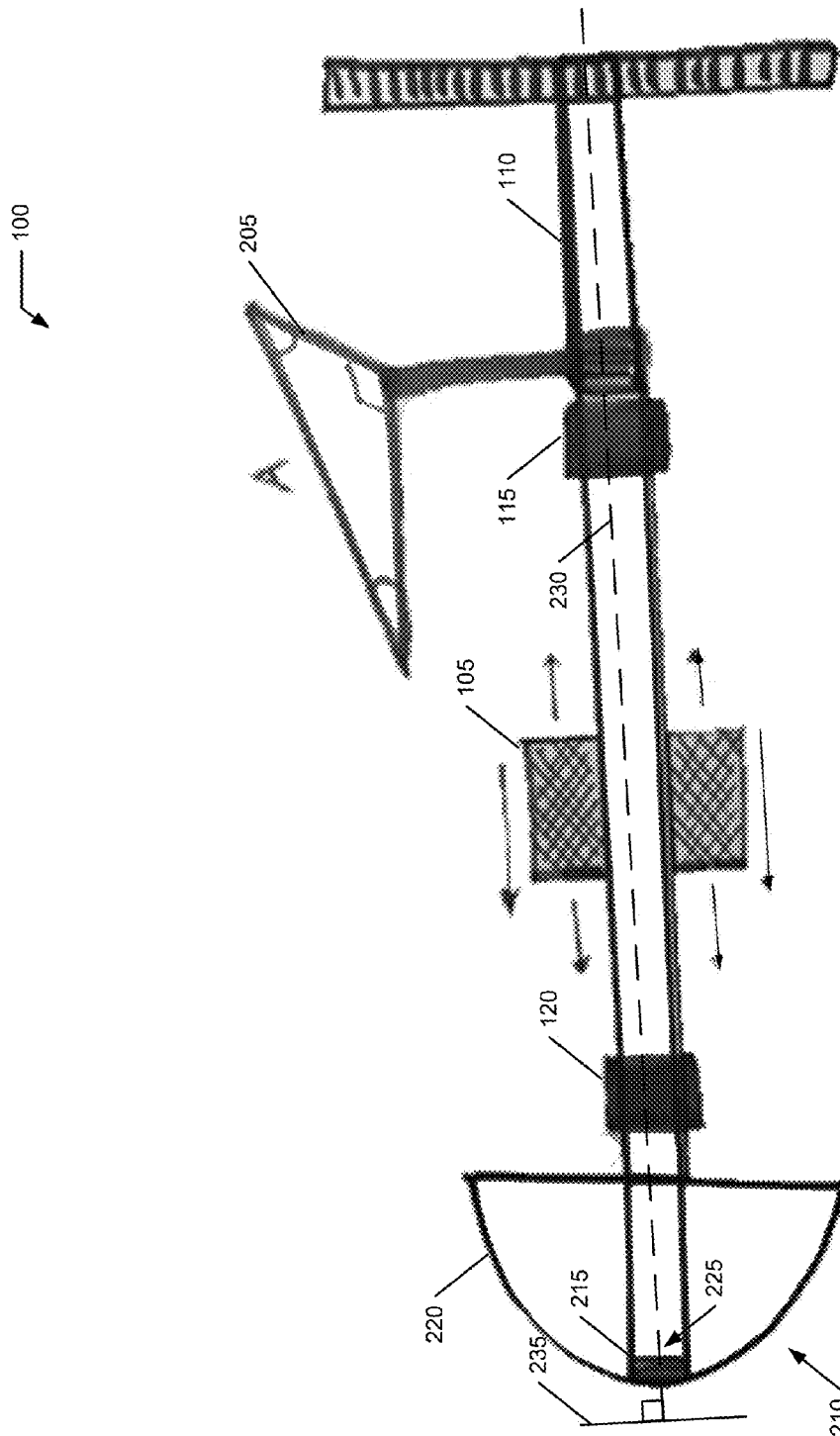

**FIG. 2**

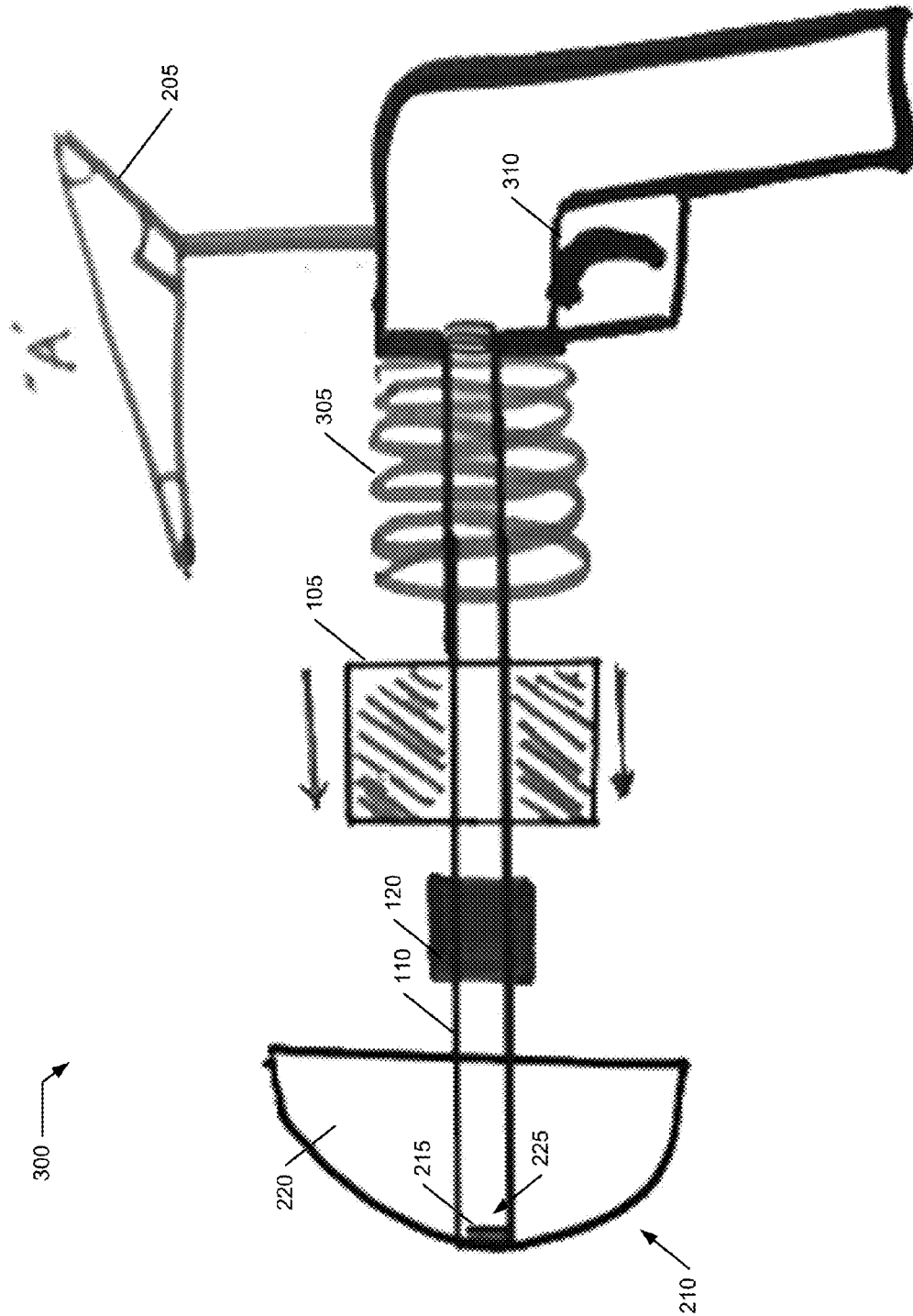


FIG. 3

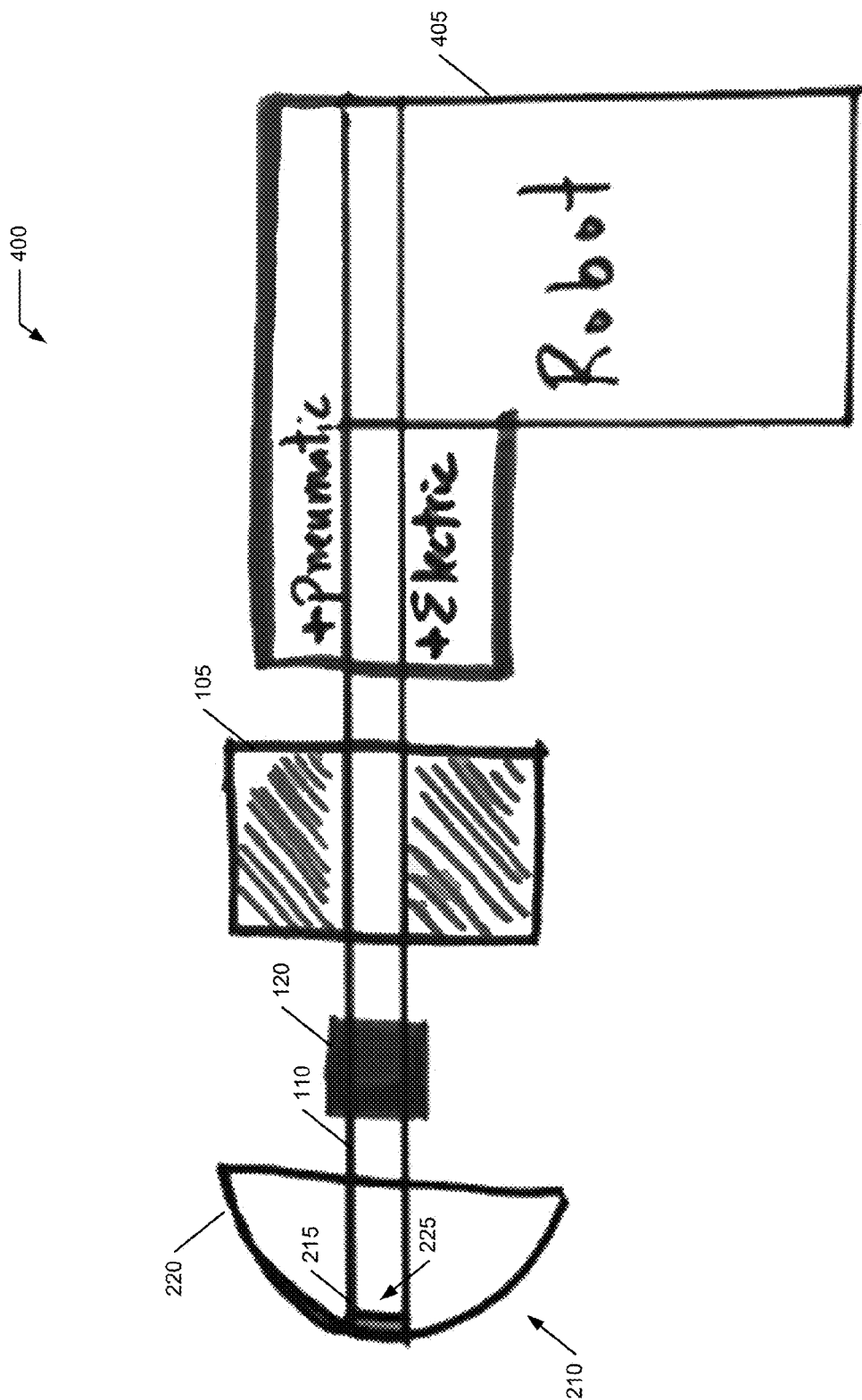


FIG. 4

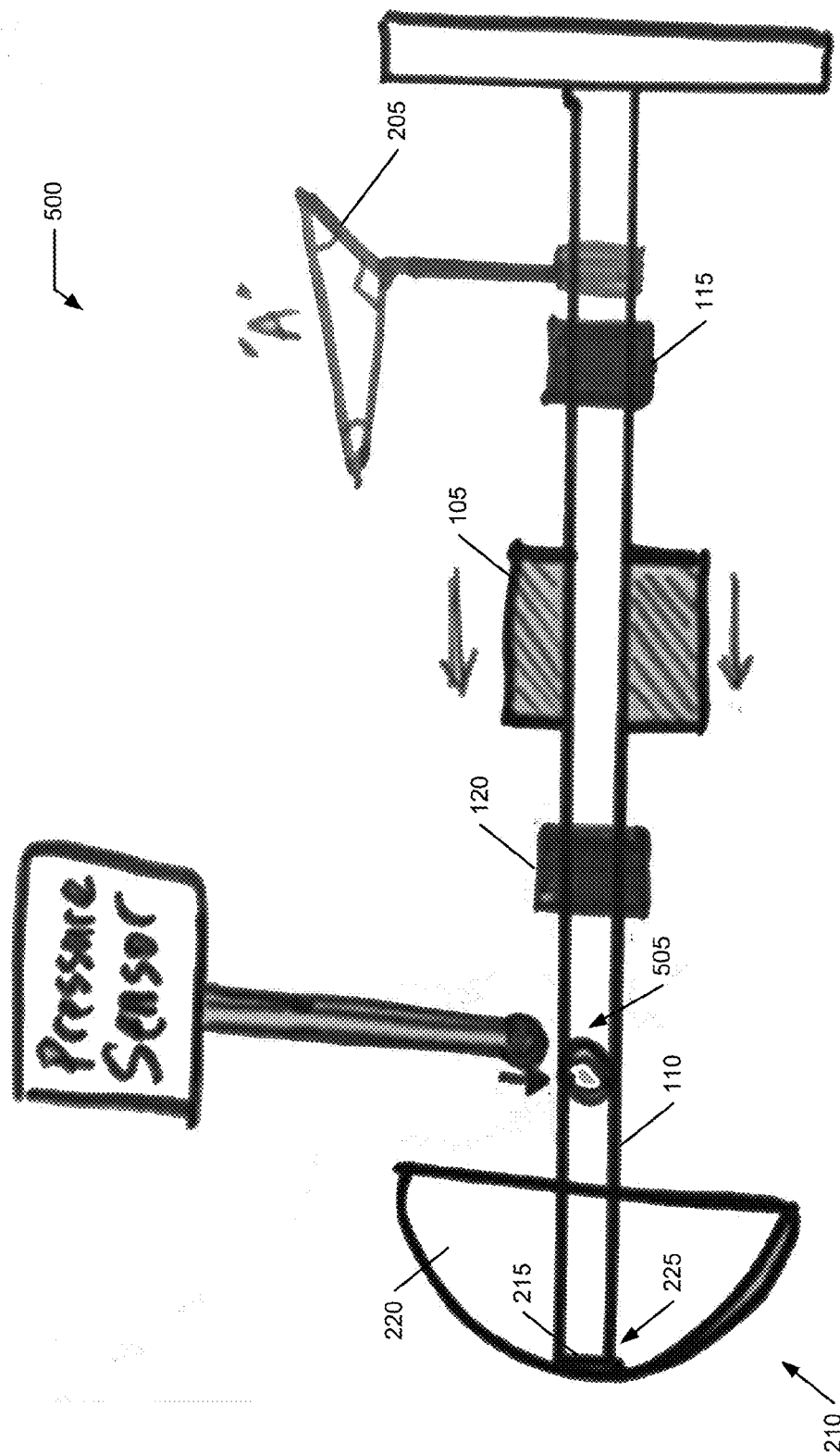
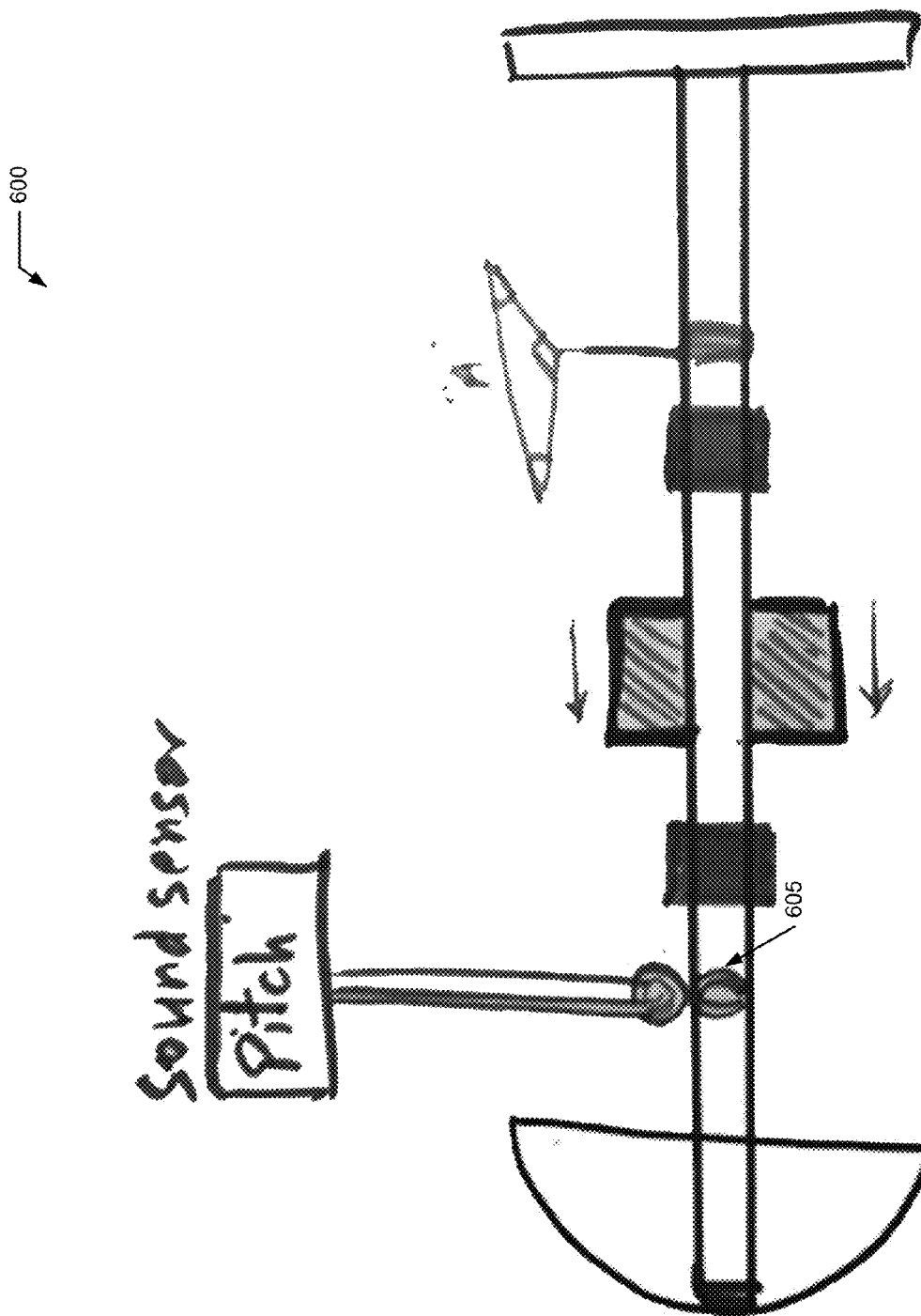


FIG. 5






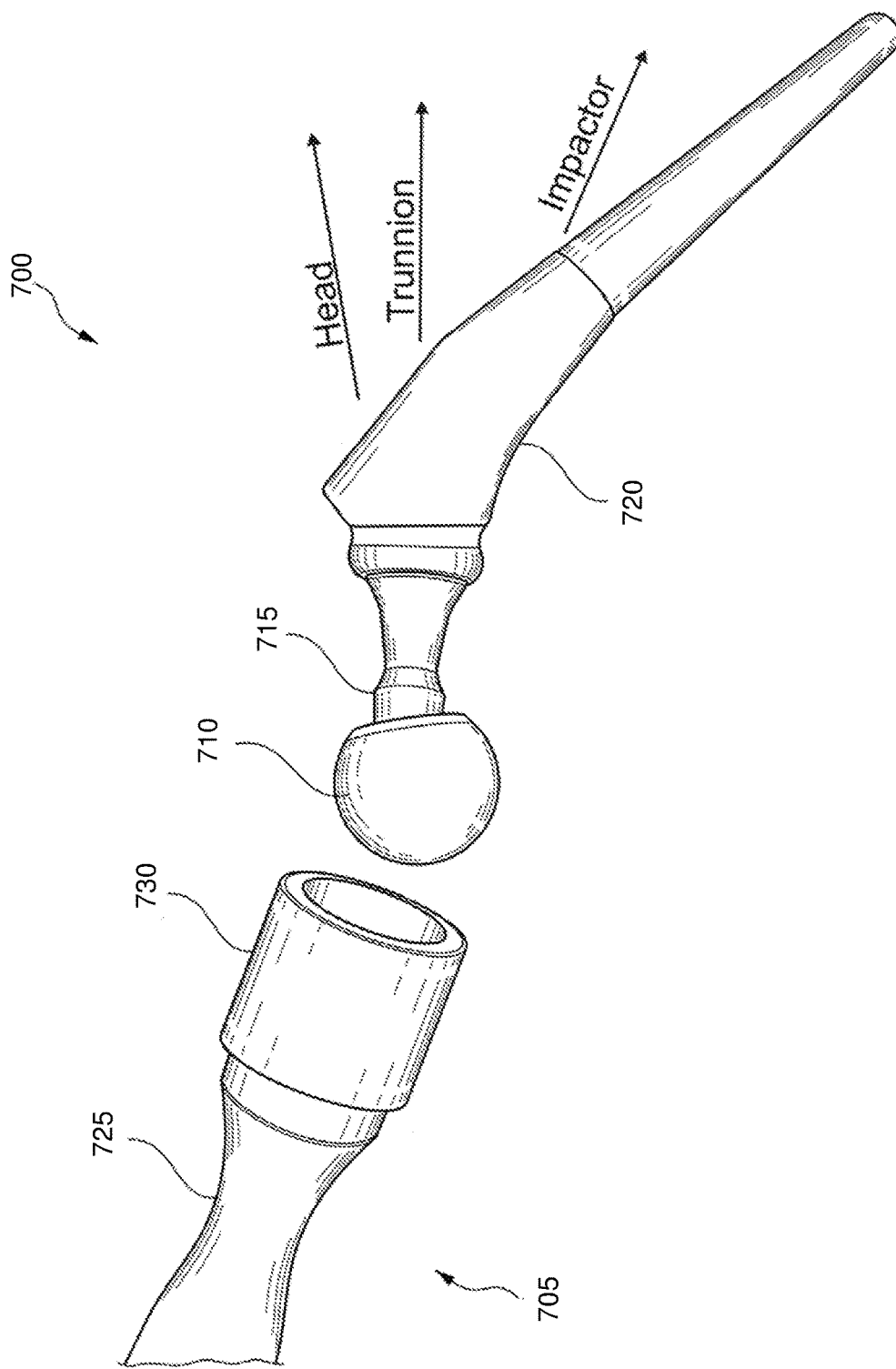


FIG. 7

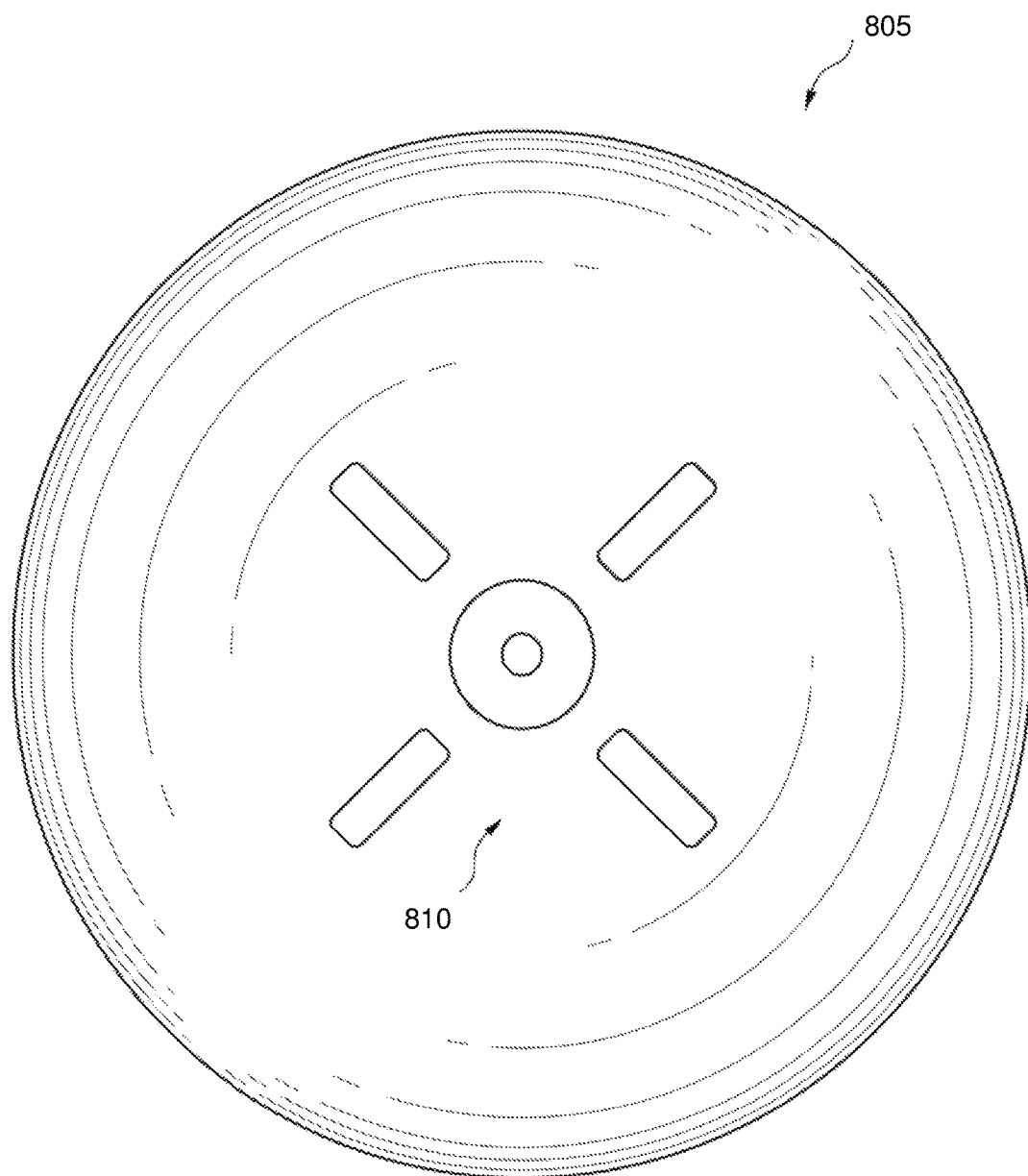

**FIG. 8**

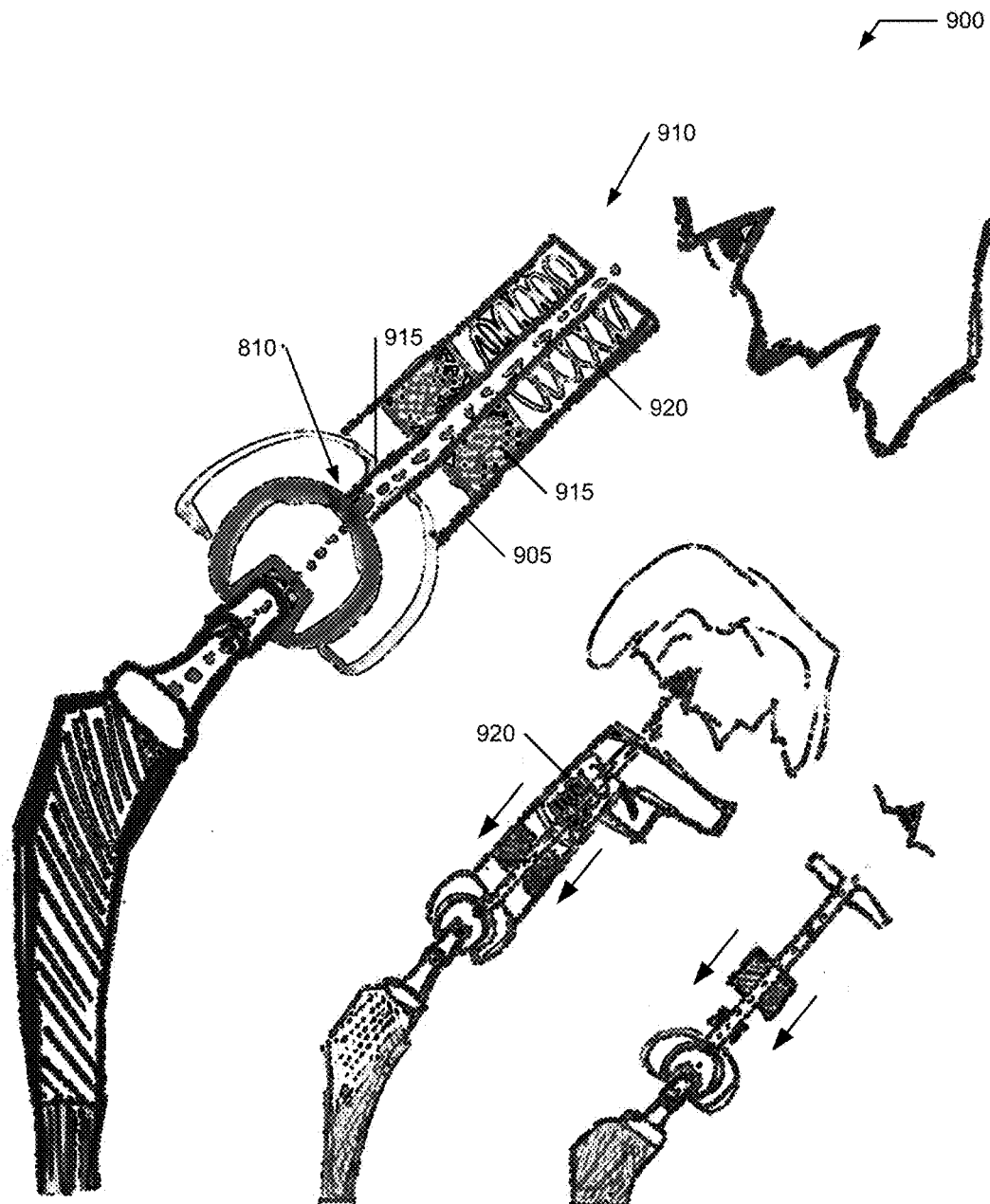


FIG. 9

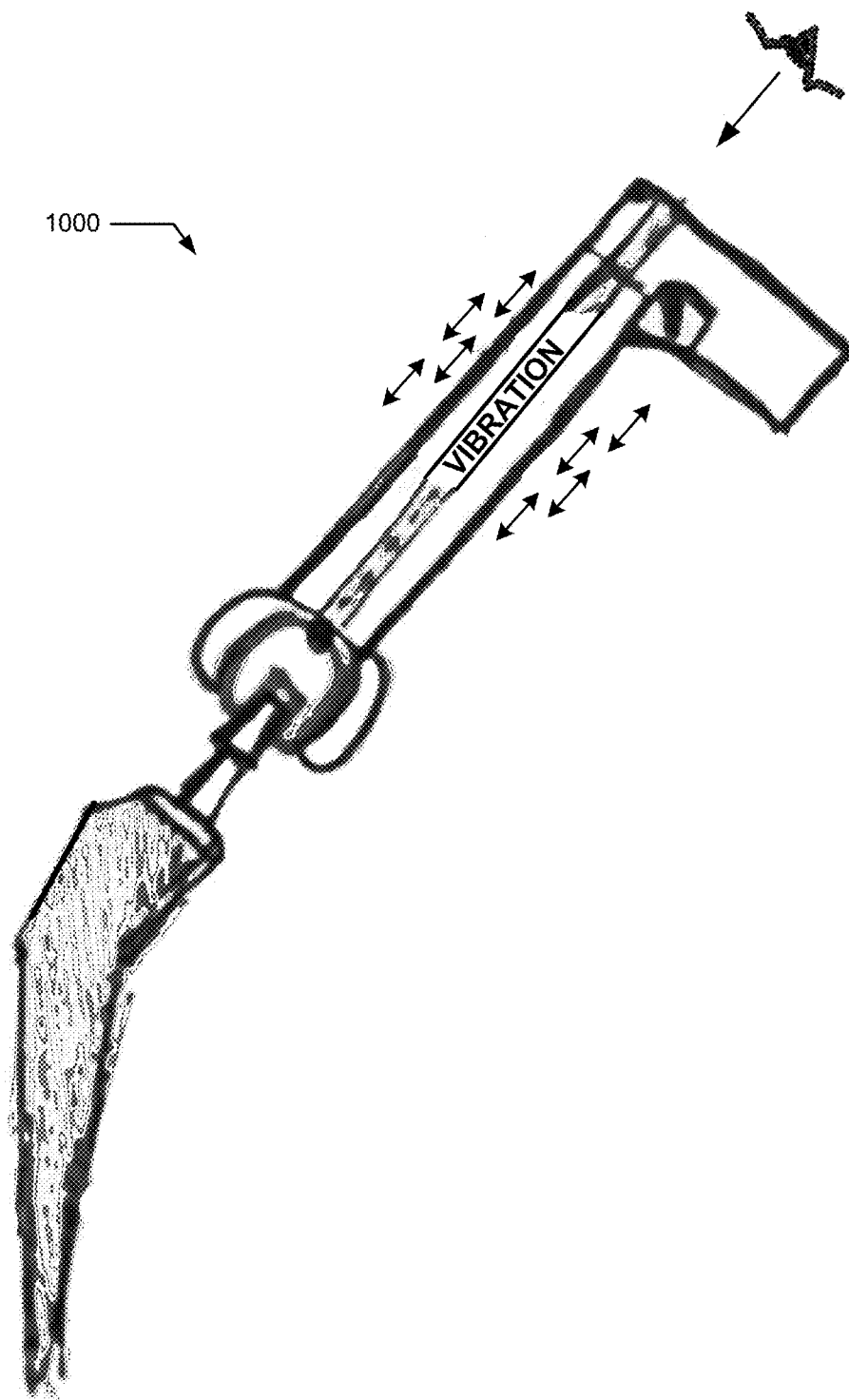


FIG. 10

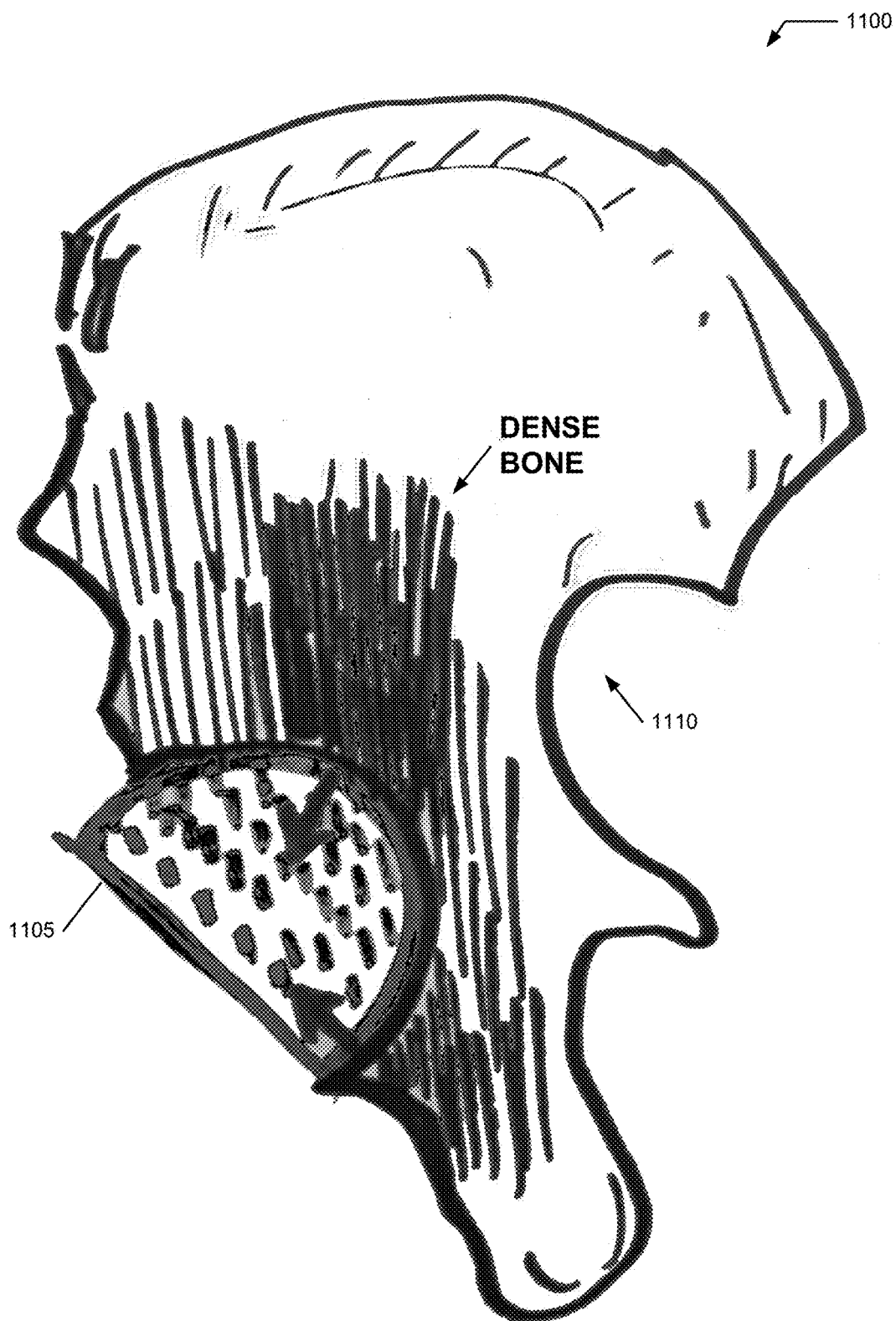
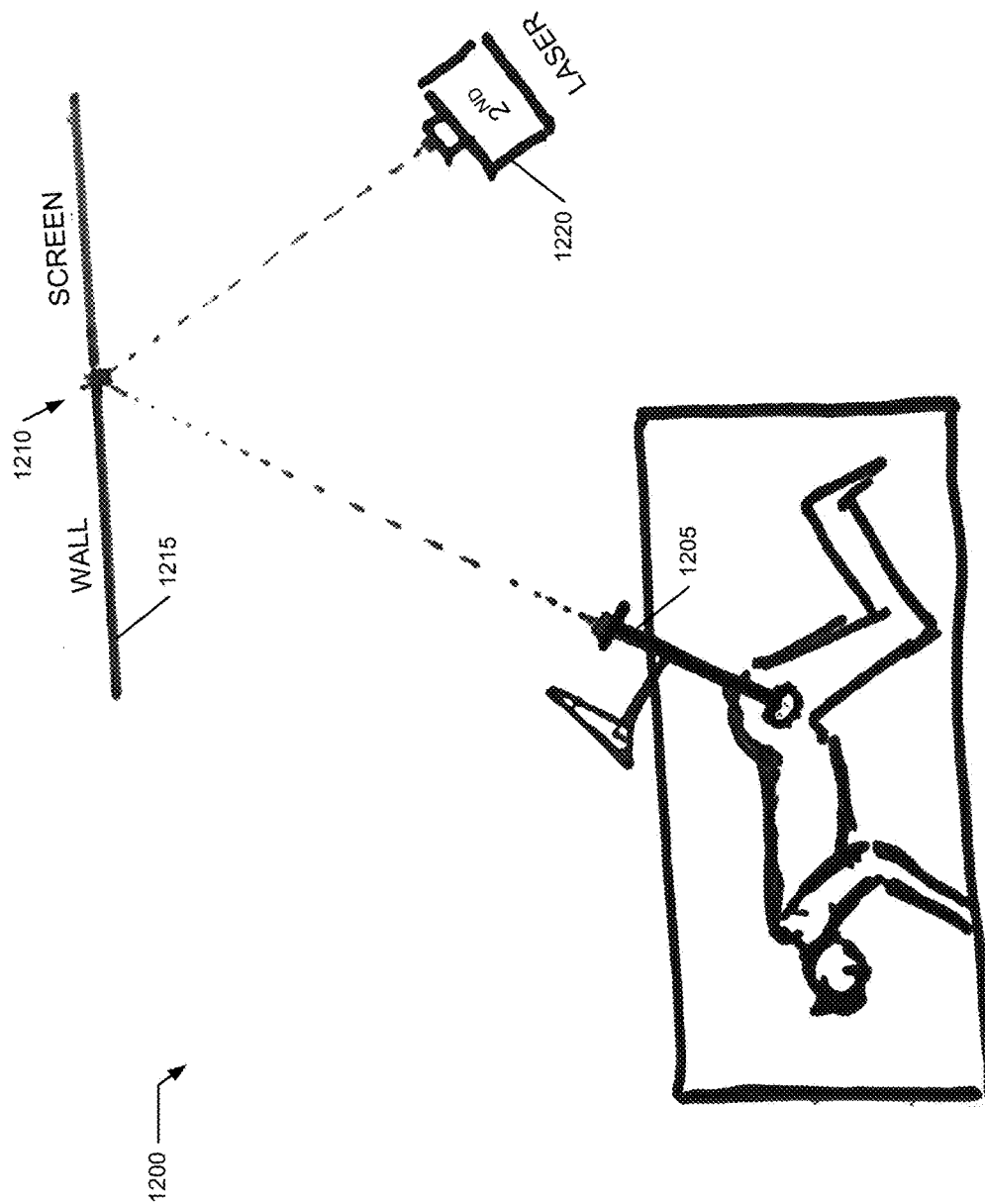


FIG. 11




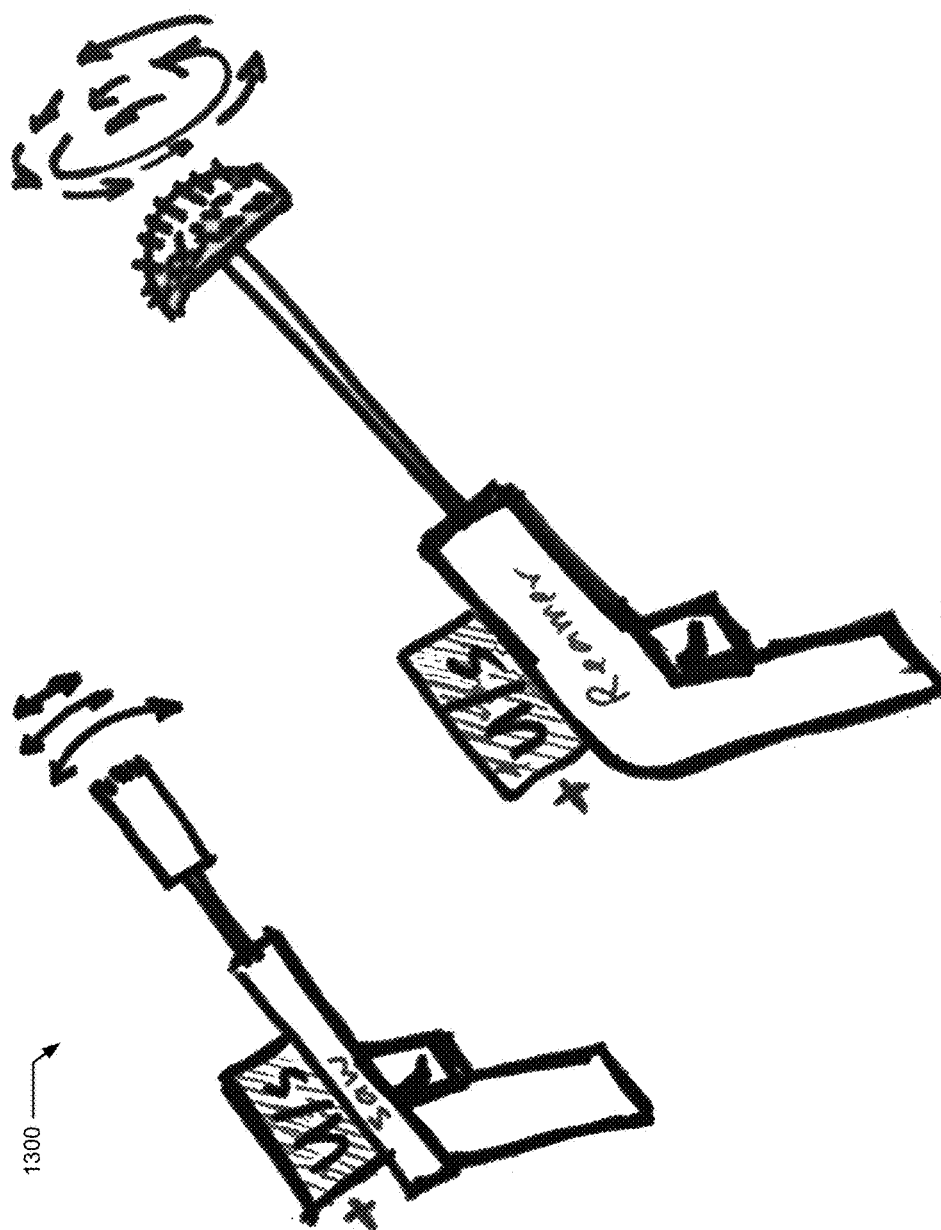



FIG. 13



**FIG. 14**

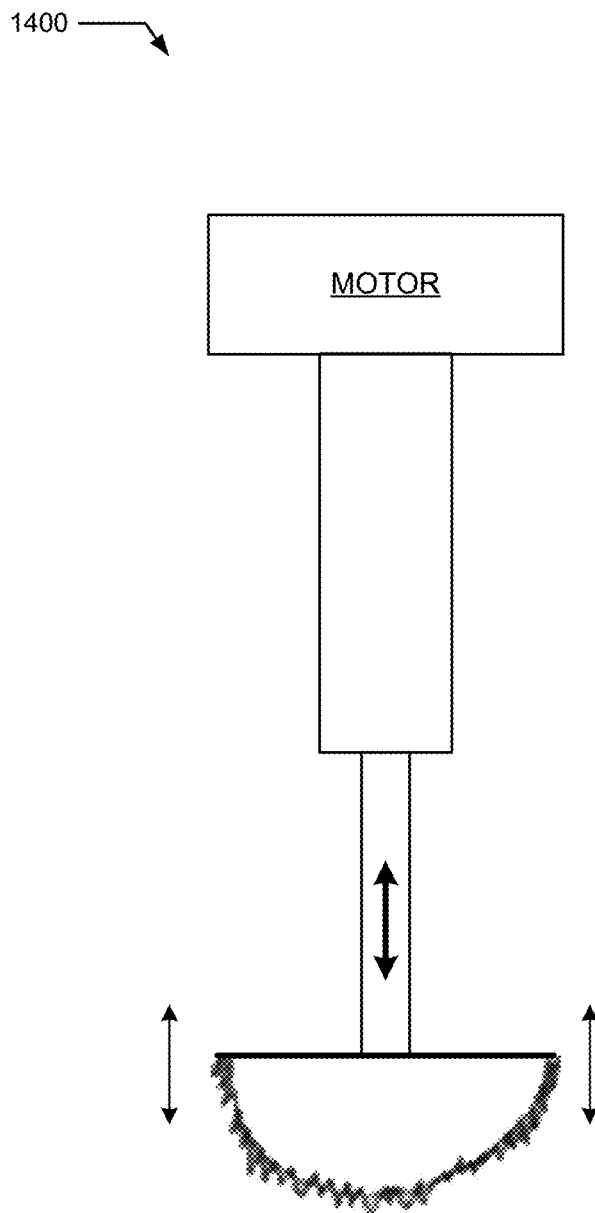

1500

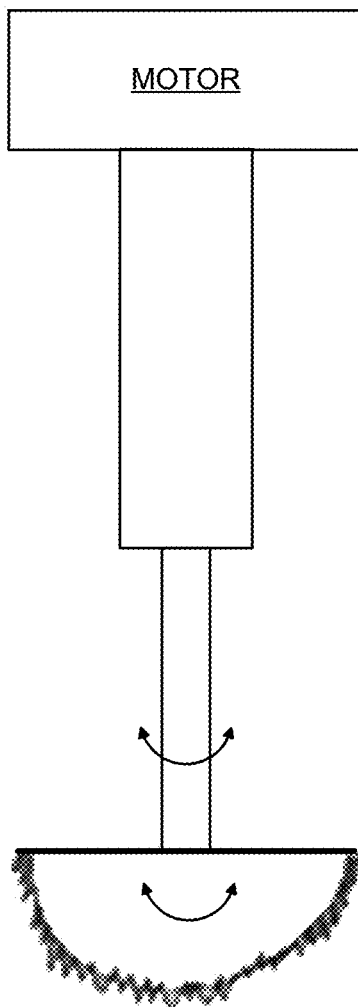

**FIG. 15**

1

**BONE PREPARATION APPARATUS AND METHOD****CROSS REFERENCE TO RELATED APPLICATIONS**

This application is a continuation-in-part of U.S. patent application Ser. No. 15/202,434 which claims benefit of U.S. Patent Application No. 62/277,294, these applications are hereby expressly incorporated by reference in their entireties for all purposes.

**FIELD OF THE INVENTION**

The present invention relates generally to installation of a prosthesis, and more specifically, but not exclusively, to improvements in prosthesis placement and positioning.

**BACKGROUND OF THE INVENTION**

The subject matter discussed in the background section should not be assumed to be prior art merely as a result of its mention in the background section. Similarly, a problem mentioned in the background section or associated with the subject matter of the background section should not be assumed to have been previously recognized in the prior art. The subject matter in the background section merely represents different approaches, which in and of themselves may also be inventions.

Earlier patents issued to the present applicant have described problems associated with prosthesis installation, for example acetabular cup placement in total hip replacement surgery. See U.S. Pat. Nos. 9,168,154 and 9,220,612, which are hereby expressly incorporated by reference thereto in their entireties for all purposes. Even though hip replacement surgery has been one of the most successful operations, it continues to be plagued with a problem of inconsistent acetabular cup placement. Cup mal-positioning is the single greatest cause of hip instability, a major factor in polyethylene wear, osteolysis, impingement, component loosening and the need for hip revision surgery.

These incorporated patents explain that the process of cup implantation with a mallet is highly unreliable and a significant cause of this inconsistency. The patents note two specific problems associated with the use of the mallet. First is the fact that the surgeon is unable to consistently hit on the center point of the impaction plate, which causes undesirable torques and moment arms, leading to mal-alignment of the cup. Second, is the fact that the amount of force utilized in this process is non-standardized.

In these patents there is presented a new apparatus and method of cup insertion which uses an oscillatory motion to insert the prosthesis. Prototypes have been developed and continue to be refined, and illustrate that vibratory force may allow insertion of the prosthesis with less force, as well, in some embodiments, of allowing simultaneous positioning and alignment of the implant.

There are other ways of breaking down of the large undesirable, torque-producing forces associated with the discrete blows of the mallet into a series of smaller, axially aligned controlled taps, which may achieve the same result incrementally, and in a stepwise fashion to those set forth in the incorporated patents, (with regard to, for example, cup insertion without unintended divergence).

There are two problems that may be considered independently, though some solutions may address both in a single solution. These problems include i) undesirable and unpre-

2

dictable torques and moment arms that are related to the primitive method currently used by surgeons, which involves manually banging the mallet on an impaction plate mated to the prosthesis and ii) non-standardized and essentially uncontrolled and unquantized amounts of force utilized in these processes.

What is needed is a system and method for improving installation of a prosthesis.

**BRIEF SUMMARY OF THE INVENTION**

Disclosed is a system and method for improving installation of a prosthesis. The following summary of the invention is provided to facilitate an understanding of some of the technical features related to prosthesis assembly and installation, and is not intended to be a full description of the present invention. A full appreciation of the various aspects of the invention can be gained by taking the entire specification, claims, drawings, and abstract as a whole. The present invention is applicable to other prosthesis in addition to acetabular cups, other modular prosthesis in addition to assembly of modular femoral and humeral prosthesis, and to other alignment and navigation systems in addition to referenced light guides.

An embodiment of the present invention may include axial alignment of force transference, such as, for example, an axially sliding hammer moving between stops to impart a non-torquing installation force. There are various ways of motivating and controlling the sliding hammer, including a magnitude of transferred force. Optional enhancements may include pressure and/or sound sensors for gauging when a desired depth of implantation has occurred.

Other embodiments include adaptation of various devices for accurate assembly of modular prostheses, such as those that include a head accurately impacted onto a trunion taper that is part of a stem or other element of the prosthesis.

Still other embodiments include an alignment system to improve site preparation, such as, for example, including a projected visual reference of a desired orientation of a tool and then having that reference marked and available for use during operation of the tool to ensure that the alignment remains proper throughout its use, such as during a reaming operation.

Further embodiments include enhancement of various tools, such as those used for cutting, trimming, drilling, and the like, with ultrasonic enhancement to make the device a better cutting, trimming, drilling, etc. device to enable its use with less strength and with improved accuracy.

Embodiments disclosed herein may include selective operational directionality. That is, for a BMD that includes vibration, it may be advantageous to control whether that vibration is driven unidirectionally and/or bidirectionally. For example, for an installation tool that installs a prosthesis into bone, it may be advantageous when a net vibratory motion is driven towards the installation site (moves toward installation) and not driven away from the installation (moves toward extraction). In a revision tool, such as disclosed in U.S. patent application Ser. No. 15/092,384, which is hereby expressly incorporated by reference in its entirety for all purposes, where it may be desired to remove a previously installed prosthesis, reversing the drive direction of the unidirectional operation helps to remove the prosthesis by providing net extractive forces on the prosthesis to be removed. As further described herein, in some implementations, it may be desirable to drive a tool operated by the BMD with a bidirectional motion. Such a system may

3

be used with a new acetabular broach, particularly with bidirectional vibratory motion.

Any of the embodiments described herein may be used alone or together with one another in any combination. Inventions encompassed within this specification may also include embodiments that are only partially mentioned or alluded to or are not mentioned or alluded to at all in this brief summary or in the abstract. Although various embodiments of the invention may have been motivated by various deficiencies with the prior art, which may be discussed or alluded to in one or more places in the specification, the embodiments of the invention do not necessarily address any of these deficiencies. In other words, different embodiments of the invention may address different deficiencies that may be discussed in the specification. Some embodiments may only partially address some deficiencies or just one deficiency that may be discussed in the specification, and some embodiments may not address any of these deficiencies.

Other features, benefits, and advantages of the present invention will be apparent upon a review of the present disclosure, including the specification, drawings, and claims.

#### BRIEF DESCRIPTION OF THE DRAWINGS

The accompanying figures, in which like reference numerals refer to identical or functionally-similar elements throughout the separate views and which are incorporated in and form a part of the specification, further illustrate the present invention and, together with the detailed description of the invention, serve to explain the principles of the present invention.

FIG. 1-FIG. 6 illustrate embodiments including installation of a prosthesis, including installation into living bone;

FIG. 1 illustrates an embodiment of the present invention for a sliding impact device;

FIG. 2 illustrates a lengthwise cross-section of the embodiment illustrated in FIG. 1 including an attachment of a navigation device;

FIG. 3 illustrates a cockup mechanical gun embodiment, an alternative embodiment to the sliding impact device illustrated in FIG. 1 and FIG. 2;

FIG. 4 illustrates an alternative embodiment to the devices of FIG. 1-3 including a robotic structure;

FIG. 5 illustrates an alternative embodiment to the devices of FIG. 1-4 including a pressure sensor to provide feedback;

FIG. 6 illustrates an alternative embodiment to the feedback system of FIG. 5 including a sound sensor to provide feedback for the embodiments of FIG. 1-5;

FIG. 7-FIG. 10 illustrate prosthesis assembly embodiments including use of variations of the prosthesis installation embodiments of FIG. 1-FIG. 6, such as may be used to reduce a risk of trunionosis;

FIG. 7 illustrates a modular prosthesis and assembly tools;

FIG. 8 illustrates a femoral head to be assembled onto a trunion attached to a femoral stem;

FIG. 9 illustrates alignment of an installation device with the femoral head for properly aligned impaction onto the trunion, such as an embodiment of FIG. 1-FIG. 6 adapted for this application;

FIG. 10 illustrates use of a modified vibratory system for assembly of the modular prosthesis;

FIG. 11-FIG. 12 illustrate an improvement to site preparation for an installation of a prosthesis;

4

FIG. 11 illustrates an environment in which a prosthesis is installed highlighting problem with site preparation; and

FIG. 12 illustrates an alignment system for preparation and installation of a prosthesis;

FIG. 13 illustrates modified surgical devices incorporating vibratory energy as at least an aid to mechanical preparation;

FIG. 14 illustrates a BMD having bidirectional longitudinal motion; and

FIG. 15 illustrates a BMD having bidirectional rotational motion.

#### DETAILED DESCRIPTION OF THE INVENTION

Embodiments of the present invention provide a system and method for improving installation of a prosthesis. The following description is presented to enable one of ordinary skill in the art to make and use the invention and is provided in the context of a patent application and its requirements.

Various modifications to the preferred embodiment and the generic principles and features described herein will be readily apparent to those skilled in the art. Thus, the present invention is not intended to be limited to the embodiment shown but is to be accorded the widest scope consistent with the principles and features described herein.

#### Definitions

Unless otherwise defined, all terms (including technical and scientific terms) used herein have the same meaning as commonly understood by one of ordinary skill in the art to which this general inventive concept belongs. It will be further understood that terms, such as those defined in commonly used dictionaries, should be interpreted as having a meaning that is consistent with their meaning in the context of the relevant art and the present disclosure, and will not be interpreted in an idealized or overly formal sense unless expressly so defined herein.

The following definitions apply to some of the aspects described with respect to some embodiments of the invention. These definitions may likewise be expanded upon herein.

As used herein, the term “or” includes “and/or” and the term “and/or” includes any and all combinations of one or more of the associated listed items. Expressions such as “at least one of,” when preceding a list of elements, modify the entire list of elements and do not modify the individual elements of the list.

As used herein, the singular terms “a,” “an,” and “the” include plural referents unless the context clearly dictates otherwise. Thus, for example, reference to an object can include multiple objects unless the context clearly dictates otherwise.

Also, as used in the description herein and throughout the claims that follow, the meaning of “in” includes “in” and “on” unless the context clearly dictates otherwise. It will be understood that when an element is referred to as being “on” another element, it can be directly on the other element or intervening elements may be present therebetween. In contrast, when an element is referred to as being “directly on” another element, there are no intervening elements present.

As used herein, the term “set” refers to a collection of one or more objects. Thus, for example, a set of objects can include a single object or multiple objects. Objects of a set also can be referred to as members of the set. Objects of a

set can be the same or different. In some instances, objects of a set can share one or more common properties.

As used herein, the term “adjacent” refers to being near or adjoining. Adjacent objects can be spaced apart from one another or can be in actual or direct contact with one another. In some instances, adjacent objects can be coupled to one another or can be formed integrally with one another.

As used herein, the terms “connect,” “connected,” and “connecting” refer to a direct attachment or link. Connected objects have no or no substantial intermediary object or set of objects, as the context indicates.

As used herein, the terms “couple,” “coupled,” and “coupling” refer to an operational connection or linking. Coupled objects can be directly connected to one another or can be indirectly connected to one another, such as via an intermediary set of objects.

The use of the term “about” applies to all numeric values, whether or not explicitly indicated. This term generally refers to a range of numbers that one of ordinary skill in the art would consider as a reasonable amount of deviation to the recited numeric values (i.e., having the equivalent function or result). For example, this term can be construed as including a deviation of  $\pm 10$  percent of the given numeric value provided such a deviation does not alter the end function or result of the value. Therefore, a value of about 1% can be construed to be a range from 0.9% to 1.1%.

As used herein, the terms “substantially” and “substantial” refer to a considerable degree or extent. When used in conjunction with an event or circumstance, the terms can refer to instances in which the event or circumstance occurs precisely as well as instances in which the event or circumstance occurs to a close approximation, such as accounting for typical tolerance levels or variability of the embodiments described herein.

As used herein, the terms “optional” and “optionally” mean that the subsequently described event or circumstance may or may not occur and that the description includes instances where the event or circumstance occurs and instances in which it does not.

As used herein, the term “size” refers to a characteristic dimension of an object. Thus, for example, a size of an object that is spherical can refer to a diameter of the object. In the case of an object that is non-spherical, a size of the non-spherical object can refer to a diameter of a corresponding spherical object, where the corresponding spherical object exhibits or has a particular set of derivable or measurable properties that are substantially the same as those of the non-spherical object. Thus, for example, a size of a non-spherical object can refer to a diameter of a corresponding spherical object that exhibits light scattering or other properties that are substantially the same as those of the non-spherical object. Alternatively, or in conjunction, a size of a non-spherical object can refer to an average of various orthogonal dimensions of the object. Thus, for example, a size of an object that is a spheroidal can refer to an average of a major axis and a minor axis of the object. When referring to a set of objects as having a particular size, it is contemplated that the objects can have a distribution of sizes around the particular size. Thus, as used herein, a size of a set of objects can refer to a typical size of a distribution of sizes, such as an average size, a median size, or a peak size.

Embodiments of the present invention may include one of more solutions to the above problems. The incorporated U.S. Pat. No. 9,168,154 includes a description of several embodiments, sometimes referred to herein as a BMD3 device, some of which illustrate a principle for breaking down large forces associated with the discrete blows of a mallet into a

series of small taps, which in turn perform similarly in a stepwise fashion while being more efficient and safer. The BMD3 device produces the same displacement of the implant without the need for the large forces from the repeated impacts from the mallet. The BMD3 device may allow modulation of force required for cup insertion based on bone density, cup geometry, and surface roughness. Further, a use of the BMD3 device may result in the acetabulum experiencing less stress and deformation and the implant may experience a significantly smoother sinking pattern into the acetabulum during installation. Some embodiments of the BMD3 device may provide a superior approach to these problems, however, described herein are two problems that can be approached separately and with more basic methods as an alternative to, or in addition to, a BMD3 device. An issue of undesirable torques and moment arms is primarily related to the primitive method currently used by surgeons, which involves manually banging the mallet on the impaction plate. The amount of force utilized in this process is also non-standardized and somewhat out of control.

With respect to the impaction plate and undesirable torques, an embodiment of the present invention may include a simple mechanical solution as an alternative to some BMD3 devices, which can be utilized by the surgeon's hand or by a robotic machine. A direction of the impact may be directed or focused by any number of standard techniques (e.g., A-frame, C-arm or navigation system). Elsewhere described herein is a refinement of this process by considering directionality in the reaming process, in contrast to only considering it just prior to impaction. First, we propose to eliminate the undesirable torques by delivering the impacts by a sledgehammer device or a structure (e.g., hollow cylindrical mass) that travels over a stainless rod.

FIG. 1 illustrates an embodiment of the present invention for a sliding impact device 100, and FIG. 2 illustrates a lengthwise cross-section of sliding impact device 100 including an attachment of a navigation device 205.

Device 100 includes a moveable hammer 105 sliding axially and freely along a rod 110. Rod 110 includes a proximal stop 115 and distal stop 120. These stops that may be integrated into rod 110 to allow transference of force to rod 110 when hammer 105 strikes distal stop 120. At a distal end 210 of rod 110, device 100 includes an attachment system 215 for a prosthesis 220. For example, when prosthesis 220 includes an acetabular cup having a threaded cavity 225, attachment system 215 may include a complementary threaded structure that screws into threaded cavity 225. The illustrated design of device 100 allows only a perfect axial force to be imparted. The surgeon cannot deliver a blow to the edge of an impaction plate. Therefore the design of this instrument is in and of itself protective, eliminating a problem of “surgeon's mallet hitting on the edge of the impaction plate” or other mis-aligned force transference, and creating undesirable torques, and hence unintentional mal-alignment of prosthesis 220 from an intended position/orientation.

A longitudinal axis 230 extends through the ends of rod 110. Attachment system 215 aligns prosthesis 220 to axis 230 when rod 110 is coupled to threaded cavity 225. An apex of prosthesis 220 (when it generally defines a hollow semi-spherical shell) supports a structure that defines threaded cavity 225 and that structure may define a plane 235 that may be tangent to the apex, with plane 235 about perpendicular to axis 230 when rod 110 engages prosthesis 220. Operation of device 100 is designed to deliver only axial (e.g., aligned with axis 230 and thus non-torqueing) forces

to prosthesis 220. Other embodiments illustrated in FIG. 3-FIG. 6 may be similarly configured.

FIG. 3 illustrates a cockup mechanical gun 300 embodiment, an alternative embodiment to the sliding impact device illustrated in FIG. 1 and FIG. 2. An alternate embodiment includes cockup mechanical gun 300 that uses the potential energy of a cocked up spring 305 to create an axially aligned impaction force. Hammer 105 is drawn back and spring 305 is locked until an operator actuates a trigger 310 to release spring 305 and drive hammer 105 along rod 110 to strike distal stop 120 and transfer an axially aligned impacting force to prosthesis 220.

Each pull of trigger 310 creates the same predetermined fixed unit of force (some alternatives may provide a variably predetermined force). The surgeon cannot deliver a misaligning impact to an impaction plate with this design.

FIG. 4 illustrates an alternative robotic device 400 embodiment to the devices of FIG. 1-3 including a robotic control structure 405. For example, device 100 and/or device 300 may be mounted with robot control structure 405 and the co-axial impacts may be delivered mechanically by a robotic tool using pneumatic or electric energy.

FIG. 5 illustrates an alternative embodiment 500 to the devices of FIG. 1-4 including a pressure sensor 505 to provide feedback during installation. With respect to management of the force required for some of these tasks, it is noted that with current techniques (the use of the mallet) the surgeon has no indication of how much force is being imparted onto the implant and/or the implant site (e.g., the pelvis). Laboratory tests may be done to estimate what range of force should be utilized in certain age groups (as a rough guide) and then fashioning a device 500, for example a modified sledgehammer 100 or cockup gun 300 to produce just the right amount of force. Typically the surgeon may use up to 2000 N to 3000 N of force to impact a cup into the acetabular cavity. Also, since some embodiments cannot deliver the force in an incremental fashion as described in association with the BMD3 device, device 500 includes a stopgap mechanism. Some embodiments of the BMD3 device have already described the application of a sensor in the body of the impaction rod. Device 500 includes sensing system/assembly 505 embedded in device 500, for example proximate rod 110 near distal end 210, and used to provide valuable feedback information to the surgeon. Pressure sensor 505 can let the surgeon know when the pressures seems to have maximized, whether used for the insertion of an acetabular cup, or any other implant including knee and shoulder implants and rods used to fix tibia and femur fractures. When pressure sensor 505 is not showing an advance or increase in pressure readings and has plateaued, the surgeon may determine it is time to stop operation/impacting. An indicator, for example an alarm can go off or a red signal can show when maximal peak forces are repeatedly achieved. As noted above, the incorporated patents describe a presence of a pressure sensor in an installation device, the presence of which was designed as part of a system to characterize an installation pulse pattern communicated by a pulse transfer assembly. The disclosure here relates to a pressure sensor provided not to characterize the installation pulse pattern but to provide an in situ feedback mechanism to the surgeon as to a status of the installation, such as to reduce a risk of fracturing the installation site. Some embodiments may also employ this pressure sensor for multiple purposes including characterization of an applied pulse pattern such as, for example, when the device includes automated control of an impacting engine coupled to the hammer. Other embodiments of this invention may

dispose the sensor or sensor reading system within a handle or housing of the device rather than in the central rod or shaft.

FIG. 6 illustrates an alternative device 600 embodiment to the feedback system of FIG. 5 including a sound sensor 605 to provide feedback for the embodiments of FIG. 1-5. Surgeons frequently use a change in pitch (sound) to gauge whether an implant (e.g., the cup) has “bottomed out” (an evaluation of a “seatedness” of the implant) and device 600 includes sound sensor 605 either attached or coupled to rod 110 or otherwise disposed separately in the operating room. Sound sensor system/assembly 605 may be used in lieu of, or in addition to, pressure sensor system/assembly 505 illustrated in FIG. 5.

FIG. 7-FIG. 10 illustrate prosthesis assembly embodiments including use of variations of the prosthesis installation embodiments of FIG. 1-FIG. 6, such as may be used to reduce a risk of trunionosis or for other advantage. FIG. 7 illustrates a modular prosthesis 700 and assembly tool 705. Prosthesis 700 includes a head 710 and a trunion taper 715 at an end of a stem 720 (e.g., a femoral stem for supporting a ball head to fit within an acetabular cup used in a total hip replacement procedure). During the procedure, the surgeon assembles prosthesis 700 by using tool 705 which may include an impact rod 725 attached to a head coupler 730. The surgeon uses tool 705 to drive head 710 onto trunion taper 715 which conventionally includes a free mallet striking tool 705. Such a procedure may be prone to the similar problems as installation of a prosthesis into an implant site, namely application of off-axis torquing forces and an uncertainty of applied force and completion of assembly.

It is believed that even a 0.1 degree mal-alignment on head 710 on trunion taper 715 may lead to progressive wear and metallosis. Variations of the embodiments of devices illustrated in FIG. 1-FIG. 6 and its associated content may be developed to help resolve this problem. In the case of “non-torquing axiallity” of forces from an assembly device, a bore of the head may define an axis, the trunion taper may define an axis, with the assembly device aligning these axes and then applying its forces in co-axial alignment with these co-axially aligned axes. Such an embodiment may reduce or eliminate any force-responsive rotations of the head with respect to the taper as the head is seated into position by the assembly device.

FIG. 8 illustrates a femoral head 805, a variation of head 710 illustrated in FIG. 7, to be assembled onto trunion taper 715 that is coupled to femoral stem 720. A center dot 810 may be placed on femoral (or humeral) head 805 to be impacted using tool 705.

FIG. 9 illustrates alignment of an installation device 900, a variation of any of devices 100-600, with femoral head 805 for properly aligned impaction onto trunion taper 715, such as an embodiment of FIG. 1-FIG. 6 adapted for this application. Such adaptation may include, for example, an axial channel 910 to view dot 810 through a slot 915, and align force transference, prior to operation of hammer 105. A sledgehammer 920 is coupled to a cock-up spring 925.

Dot 810 can be aligned with an impactor/device/gun. Once axial alignment, such as through the sight channel, has been confirmed, a sledgehammer, a cockup gun, or other similar device can bang the impactor onto femoral (humeral) head 805 to impact it on trunion taper 715. The co-axiality of the head and the device can be confirmed visually (for example, through a hollow cylinder that comprises a center shaft of the device) or with a variety of electronic and laser methods.

FIG. 10 illustrates use of a modified vibratory system 1000, a variation of installation device 900 for assembly of the modular prosthesis illustrated in FIG. 7. Alternatively to device 900, a variation of the BMD3 device can be used to insert the femoral and humeral heads 710 onto trunion taper 715. For example, a version of the BMD3 device where femoral head 710 is grasped by a “vibrating gun” and introduced methodically and incrementally onto trunion taper 715. Since there are no large forces being applied to the head/trunion junction, there is essentially no possibility, or a reduced possibility, of head 710 seating onto trunion taper 715 in a misaligned fashion. It would be possible to use the same technique of marking the center of head 710 and lining it up with trunion taper 715 and device axially before operating the device.

FIG. 11-FIG. 12 illustrate an improvement to site 1100 preparation for an installation of a prosthesis 1105. FIG. 11 illustrates an environment 1100 in which prosthesis 1105 is installed highlighting a problem with site preparation for a prosthesis installation procedure having variable density bone (line thickness/separation distance reflecting variable bone density) of acetabulum 1110.

There is a secondary problem with the process of acetabular preparation and implantation that leads to cup malalignment. Currently, during the process of acetabular reaming, surgeons make several assumptions. One common assumption is that the reamer is fully seated in a cavity and surrounded on all sides by bone. Another common assumption is that the bone that is being reamed is uniform in density. Imagine a carpenter that is preparing to cut a piece of wood with a saw. Now imagine that parts of this piece of wood are embedded with cement and some parts of the piece of wood are hollow and filled with air. The carpenter’s saw will not produce a precise cut on this object. Some parts are easy to cut and some parts are harder to cut. The saw blades skives and bends in undesirable ways. A similar phenomenon happens in acetabular preparation with a reamer and when performing the cuts for knee replacement with a saw. With respect to the acetabulum, the side of the cavity that is incomplete (side of the reamer that is uncovered) will offer less resistance to the reamer and therefor the reamer preferentially reams towards the direction of the uncovering. Second, the reamer cuts the soft bone much more easily than the dense and sclerotic bone, so the reamer moves away from the sclerotic bone and moves toward the soft bone. From a machining perspective, the reaming and preparation of the acetabulum may not be concentric or precise. This maybe a significant factor in the surgeon’s inability to impact the cup in the desired location.

FIG. 12 illustrates an alignment system 1200 for preparation and installation of a prosthesis to help address/minimize this effect. A first step that can be taken is to include directionality into the process of reaming at the outset, and not just at the last step during impaction. Current technique allows the surgeon to ream the cup haphazardly moving the reamer handle in all directions, being ignorantly unaware that he is actually creating a preference for the sinking path of the acetabular implant. Ultimately the direction in which the surgeon reams may in fact be determining the position/path of the final implant. The surgeon then impacts the cup using the traditional A-frame or any of the currently used intra-operative measurement techniques such as navigation or fluoroscopy. These methods provide information about the position of the cup either as it is being implanted or after the implantation has occurred. None of these techniques predetermine the cup’s path or function to guide the cup in the correct path.

Proposed is a method and a technique to eliminate/reduce this problem. Before the surgeon begins to ream the acetabulum, the reamer handle should be held, with an A-frame attached, in such a way to contemplate the final position of the reamer and hence the implant, (e.g., hold the reamer in 40 degree abduction and 20 degree anteversion reaming is started). This step could also be accomplished with navigation or fluoroscopy. The surgeon could, for example, immediately mark this position on a screen or the wall in the operating room as described below and as illustrated in FIG. 12. After the anticipated position of the reamer is marked, the surgeon can do whatever aspect of reaming that needs to be done. For example the first reaming usually requires medialization in which the reamer is directed quite vertically to ream in to the pulvinar. Typically three or four reamings are done. First, the acetabular cavity is medialized. The other reamings function to get to the subchondral bone in the periphery of the acetabulum. One solution may be that after each reaming, the reamer handle be held in the final anticipated position of the implant. In some cases it may be difficult to have an A-frame attached to every reamer and to estimate the same position of the reamer in the operating space accurately with the A-frame.

An alternative to that is also proposed to address this process. For example, at a proximal end of the reamer shaft handle will be placed a first reference system 1205, for example a laser pointer. This laser pointer 1205 will project a spot 1210 either on a wall or on a screen 1215, a known distance from the operating room table. That spot 1210 on wall 1215 (or on the screen) is then marked with another reference system 1220, for example a second independent laser pointer that sits on a steady stand in the operating room. Thereafter manipulating the shaft handle so that the first reference system has the desired relationship, example co-aligned, with the second reference system, the surgeon knows that the device attached to the handle has the desired orientation. So when the first reamer is held in the anticipated and desired final alignment of the implant (e.g., 40 degree abduction, 20 degree anteversion for many preferred installation angles of an acetabular cup), the laser pointer at the proximal end of the reamer handle projects a spot on the wall or screen. That spot is marked with the second stationary laser, and held for the duration of the case. All subsequent reamings will therefore not require an A-frame to get a sense of the proper alignment and direction of the reamer. The surgeon assures that no matter how he moves the reamer handle in the process of reaming of the acetabulum, that the reaming finishes with the reamer handle (laser pointer) pointing to the spot on the wall/screen. In this manner, directionality is assured during the reaming process. In this way the sinking path of the actual implant is somewhat predetermined. And no matter what final intra-operative monitoring technique is used (A-frame, C-Arm, Navigation) that the cup will likely seat/sink more closely to the desired final position.

FIG. 13 illustrates modified surgical devices 1300 incorporating vibratory energy as at least an aid to mechanical preparation. Also proposed herein is another concept to address a problem associated with non-concentric reaming of the acetabulum caused by variable densities of the bone and the uncovering of the reamer. Imagine the same carpenter has to cut through a construct that is made out of wood, air, and cement. The carpenter does not know anything about the variable densities of this construct. There are two different saws available: one that cuts effectively through wood only, and ineffectively through the cement. Also available is a second saw that cuts just as effectively

11

through cement as wood. Which of these saws would improve a chance of producing a more precise cut? Proposed is a mixing of ultrasonic energy with the standard oscillating saw and the standard reamer. In effect any oscillating equipment used in orthopedics, including the saw, reamer, drill, and the like may be made more precise in its ability to cut and prepare bone with the addition of ultrasonic energy. This may feel dangerous and counterintuitive to some, however, the surgeon typically applies a moderate amount of manual pressure to the saw and reamers, without being aware, which occasionally causes tremendous skiving, bending and eccentric reaming. An instrument that does not require the surgeon's manual force maybe significantly safer and as well as more precise and effective.

A further option includes disposition of a sensor in the shaft of the ultrasonic reamers and saws so that the surgeon can ascertain when hard versus soft bone is being cut, adding a measure of safety by providing a visual numerical feedback as to the amount of pressure being utilized. This improvement (the ability to cut hard and soft bone with equal efficacy) will have tremendous implications in orthopedic surgery. When the acetabular cavity is prepared more precisely, with significantly lower tolerances, especially when directionality is observed, the acetabular implant (cup) may more easily follow the intended sinking path.

Other applications of this concept could be very useful. Pressfit and ingrowth fixation in total knee replacements in particular (as well as ankle, shoulder and other joints to a lesser degree) are fraught with problems, particularly that of inconsistent bony ingrowth and fixation. The fact that a surgeon is unable to obtain precise cuts on the bone may be a significant factor in why the bone ingrowth technology has not gotten off the ground in joints other than the hip. The problem is typically blamed on the surgeon and his less than perfect hands. The experienced surgeon boasts that only he should be doing this operation (i.e.: non-cemented total knee replacement). This concept (a more precise saw that cuts hard and soft bone equally allowing lower tolerances) has huge potential in orthopedics, in that it can lead to elimination of the use of cement in orthopedic surgery altogether. This can spark off the growth and use of bone ingrowth technology in all aspects of joint replacement surgery which can lead to tremendous time saving in the operating room and better results for the patients.

FIG. 14 illustrates a BMD 1400 having bidirectional longitudinal motion; and FIG. 15 illustrates a BMD 1500 having bidirectional rotational motion. In previous discussions of BMD3 vibratory and operational devices, specific directionality controls of the movement were not addressed as described herein. Many vibratory systems are "driven" in one-direction based upon a particular application. Disclosed herein are devices that have intentionally designed and allow for, based upon application, for unidirectionality in applied force by an oscillatory engine. For procedures and processes relating to preparing an installation site, installing a prosthesis, and revising/removing an installed prosthesis, there may be advantages in different directionalities in different contexts. Rather than having three different tools, the present disclosure contemplates a tool having multiple selectable directionalities allowing it to be used in different procedures.

Also disclosed is a new type of cavity formation tool (for hip replacement in preparation of the pelvic bone) that may advantageously employ bidirectional vibratory motion: a broach for the acetabulum cavity preparation.

BMD3 bidirectional vibratory tool: The BMD3 vibratory tool was initially created and envisioned for vibratory insertion of prosthesis into bone. During the experimentation of

12

BMD3 vibratory tool we discovered that vibratory energy can be unidirectional in forward and backward directions or it can be bidirectional. We have described the effectiveness and use of unidirectional forward vibrating BMD3 tool for insertion of a prosthesis (in particular acetabular prosthesis) into bone. We now propose use of bidirectional BMD3 vibratory tool for the purpose of preparing bone, and in particular the acetabular cavity.

BMD3 bidirectional vibratory tool for preparation of bone, and in particular the acetabular cavity: The use of a Acetabular Broach: a new idea. BMD3 bi-directional vibratory tool can be used for preparation of bone (any cavity of bone that needs to be prepared for application of a prosthesis, but especially the acetabulum, as well as the proximal femur, proximal tibia, proximal humerus, and any other long bone in the body that receives a prosthesis). With regards to the acetabulum, unlike the other bones discussed above, this structure has never before been prepared with a broach, but rather always prepared with a hemispherical "cheese grater type" reamers that rotates in one direction (forward). We are proposing that the acetabulum be prepared with a broach using one of the two degrees of freedom for oscillation

(1. Longitudinal and 2. rotational), utilizing a bidirectional BMD vibratory tool. The outer surface of this broach will very closely resemble the rough surface of the prosthesis, with high coefficient of static friction. We have seen this method in action in our experiments, particularly at higher frequencies of around 300 hertz, and believe that this method of acetabular preparation will provide a cut surface that is much more precise and conferring the ability to produce lower tolerances. This method may also allow preparation of acetabular cavity in "half" sizes. Currently the cavity is reamed in 1 mm intervals. It may be much easier to prepare the acetabulum with ½ mm interval broaches than ½ mm reamers. Half size broaching may dramatically improve the ability of the surgeon to cut and prepare the acetabular precisely and at lower tolerances.

For purposes of review we recall the equation  $FR=K \cdot x \cdot Us$ . Where  $x$  is represents the amount of under reaming and the shape of the cup being inserted.

$X$  is controlled by the amount of under or over reaming of the acetabulum. In the past when the surfaces of the cup were not as rough (lower coefficient of static friction, i.e. Zimmer Fiber Metal cup), surgeons used to under ream by 2 mm. Now most companies recommend under reaming by 1 mm, since the surfaces of most cups are much more rough with better porosity characteristics that allow better and quicker bony ingrowth. Sometimes when the surgeon has difficulty seating the cup, he/she reams line to line, and describes this action as "touching up the rim". This action however, many times, eliminates the compressive quality of the acetabulum by decreasing the value of  $x$  towards zero. This issue brings attention to the problem that we have described which is that the surgeon does not have anything but a most basic understanding of the spring like qualities of bone. If he/she is can understand the basic science involved in this system, he can then use the proper tools to appropriately fine tune the pelvis for a good press fit fixation, without fear of under seating or fracture. There is a huge market need for better tools to prepare (fine tune) the acetabulum, for good press fit fixation.

Current techniques utilize 'cheese grater type' hemispherical reamers to prepare the bed of the acetabulum. As discussed in our BMD4 paper the quality of acetabular bone can be drastically different between patients and even within the same patient, particularly at different locations around the acetabular fossa. Some parts of the bone are soft, and

13

some are hard. Current cheese grater hemispherical reamers come in 1 mm intervals. This creates two specific problems: 1. The current acetabular reamers in 1 mm intervals for preparation of the acetabular bone do not provide the ability to precisely machine the acetabulum, and obtain lower tolerances, and therefore proper tuning of the pelvic bone. 2. No method exists to cut hard and soft bone with the same level of effectiveness, i.e.: hard bone always pushes the reamers towards the soft bone which ends up being chewed up more, and in that sense, a perfect hemisphere is not created with current cheese grater reaming techniques. We therefore are proposing two distinct and separate solutions which we believe can remedy this problem of poor quality acetabular preparation.

1. The creation of half reamers. The production and use of half reamers gives the surgeon the ability to ream up or down by half millimeters. Which gives him/her the ability to fine tune x more precisely, and therefore FR more precisely. This basically gives the surgeon a better set of tuning forks to obtain better tension for the acetabulum and utilize its viscoelastic properties to his/her advantage to obtain a better press fit fixation.

2. Ultrasonic assisted reaming or broaching: Lastly, we believe that there is some room for creating a better cutting tool by adding ultrasonic energy to either the acetabular broach described above or the acetabular half reamers described above to create an ultrasonic assisted reaming or broaching of the acetabulum for obtaining a more precise cut and at a lower tolerance. We believe this is a new and novel idea that can be considered for preparation of the acetabulum for obtaining better tension of the pelvis for application of an acetabular prosthesis.

The system and methods above has been described in general terms as an aid to understanding details of preferred embodiments of the present invention. In the description herein, numerous specific details are provided, such as examples of components and/or methods, to provide a thorough understanding of embodiments of the present invention. Some features and benefits of the present invention are realized in such modes and are not required in every case. One skilled in the relevant art will recognize, however, that an embodiment of the invention can be practiced without one or more of the specific details, or with other apparatus, systems, assemblies, methods, components, materials, parts, and/or the like. In other instances, well-known structures, materials, or operations are not specifically shown or described in detail to avoid obscuring aspects of embodiments of the present invention.

Reference throughout this specification to "one embodiment", "an embodiment", or "a specific embodiment" means that a particular feature, structure, or characteristic described in connection with the embodiment is included in at least one embodiment of the present invention and not necessarily in all embodiments. Thus, respective appearances of the phrases "in one embodiment", "in an embodiment", or "in a specific embodiment" in various places throughout this specification are not necessarily referring to the same embodiment. Furthermore, the particular features, structures, or characteristics of any specific embodiment of the present invention may be combined in any suitable manner with one or more other embodiments. It is to be understood that other variations and modifications of the embodiments of the present invention described and illustrated herein are possible in light of the teachings herein and are to be considered as part of the spirit and scope of the present invention.

14

It will also be appreciated that one or more of the elements depicted in the drawings/figures can also be implemented in a more separated or integrated manner, or even removed or rendered as inoperable in certain cases, as is useful in accordance with a particular application.

Additionally, any signal arrows in the drawings/Figures should be considered only as exemplary, and not limiting, unless otherwise specifically noted. Combinations of components or steps will also be considered as being noted, where terminology is foreseen as rendering the ability to separate or combine is unclear.

The foregoing description of illustrated embodiments of the present invention, including what is described in the Abstract, is not intended to be exhaustive or to limit the invention to the precise forms disclosed herein. While specific embodiments of, and examples for, the invention are described herein for illustrative purposes only, various equivalent modifications are possible within the spirit and scope of the present invention, as those skilled in the relevant art will recognize and appreciate. As indicated, these modifications may be made to the present invention in light of the foregoing description of illustrated embodiments of the present invention and are to be included within the spirit and scope of the present invention.

Thus, while the present invention has been described herein with reference to particular embodiments thereof, a latitude of modification, various changes and substitutions are intended in the foregoing disclosures, and it will be appreciated that in some instances some features of embodiments of the invention will be employed without a corresponding use of other features without departing from the scope and spirit of the invention as set forth. Therefore, many modifications may be made to adapt a particular situation or material to the essential scope and spirit of the present invention. It is intended that the invention not be limited to the particular terms used in following claims and/or to the particular embodiment disclosed as the best mode contemplated for carrying out this invention, but that the invention will include any and all embodiments and equivalents falling within the scope of the appended claims. Thus, the scope of the invention is to be determined solely by the appended claims.

What is claimed as new and desired to be protected by Letters Patent of the United States is:

1. A device for imparting a vibratory bone preparation force, comprising:

a rod having a shaft including a proximal end, a distal end spaced apart from said distal end, and a longitudinal axis extending from said proximal end to said distal end through said rod;

a motor producing a periodic motion;

a driver system, coupled to said proximal end and to said motor, producing a driven vibratory rod motion for said distal end from said periodic motion;

an attachment system coupled to said distal end, said attachment system configured to engage an attachment structure; and

a set of attachments, each particular attachment from said set of attachments including said attachment structure;

wherein said driven vibratory rod motion includes a selectable vibratory mode chosen from a vibratory mode group including a bidirectionally driven vibratory mode having a first driven vibratory direction and a second driven vibratory direction.

2. The device of claim 1 wherein said vibratory mode group further includes a first unidirectional vibratory mode and further comprising a vibratory mode selector to select a



15

particular one vibratory mode from said bidirectionally driven vibratory mode and said first unidirectional vibratory mode.

3. The device of claim 1 wherein said bidirectionally driven vibratory mode includes a bidirectional longitudinal motion vibratory mode.

4. The device of claim 3 wherein said vibratory mode group further includes a first unidirectional longitudinal vibratory mode and further comprising a vibratory mode selector to select a particular one vibratory mode from said bidirectional longitudinal motion vibratory mode and said first unidirectional longitudinal vibratory mode.

5. The device of claim 1 wherein said bidirectionally driven vibratory mode includes a bidirectional rotational motion vibratory mode.

6. The device of claim 5 wherein said vibratory mode group further includes a first unidirectional rotational vibratory mode and further comprising a vibratory mode selector to select a particular one vibratory mode from said bidirectional rotational vibratory mode and said first unidirectional rotational vibratory mode.

7. The device of claim 1 wherein said set of attachments includes a bone broach.

16

8. The device of claim 1 wherein said set of attachments includes a set of broaches including a range of varying sizes and wherein said range of varying sizes include ½ millimeter variations.

9. The device of claim 5 wherein said set of attachments includes a bone broach.

10. The device of claim 9 wherein said set of attachments includes a set of broaches including a range of varying sizes and wherein said range of varying sizes include ½ millimeter variations.

11. The device of claim 1 wherein said set of attachments includes an acetabular broach.

12. The device of claim 3 wherein said set of attachments includes an acetabular broach.

13. The device of claim 1 wherein said bidirectionally driven vibratory mode operates at an ultrasonic frequency.

14. The device of claim 9 wherein said bidirectionally driven vibratory mode operates at an ultrasonic frequency.

15. The device of claim 11 wherein said bidirectionally driven vibratory mode operates at an ultrasonic frequency.

\* \* \* \* \*